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CONTENTS

- 01 GENERAL INFORMATION
- 02 PACKAGE CONTENTS
- 03 TECHNICAL SPECIFICATIONS
- 04 DESCRIPTION OF DESIGN AND CONTROLS
- 05 USAGE
- 06 SAFETY PRECAUTIONS
- 07 POSSIBLE ISSUES AND SOLUTIONS
- 08 SYMBOLS

Before using the massager, please read the user manual carefully, paying special attention to the safety instructions and contraindications.

**Keep this manual for
future reference.**

GENERAL INFORMATION

EN

The Zabobon massager provides a three-dimensional massage with adjustable depth of impact, equipped with a six-element roller system that mimics the movements of a masseur's hands.

Shiatsu Roller Effect

The massage system replicates the popular Japanese **Shiatsu** technique, designed for deep muscle stimulation in the neck, shoulders, and upper back. The rollers perform combined movements with varying pressure intensity, delivering precise and effective action.

Infrared Heating

The built-in infrared heat source improves blood circulation, relieves pain, and promotes muscle relaxation. The heating function is manually activated.

Safe Operation

The device **automatically** shuts off after 10 minutes of continuous operation to prevent overheating. This protection is implemented at the hardware level and does not require user intervention.

Convenient Controls

All control elements are located on the front part of the massager's handle, providing **quick access** to mode switching and intensity adjustments. Factory settings are configured to optimal levels.

Cordless Operation

The rechargeable lithium-ion battery allows the massager to be used without a power outlet. Its compact design and lightweight construction make it easy to transport and suitable for use on various body areas.

Comfortable Materials

The massager's body is covered with soft and durable **eco-leather**. The working area with rollers is lined with breathable, hypoallergenic fabric that is safe for sensitive skin.

TECHNICAL SPECIFICATIONS

Power Supply

USB Type-C, 5 V, 2 A

Power:

10 W

Battery

**Lithium-ion (Li-Ion), 7.4 V,
2000 mAh (14.8 Wh)**

Heating Temperature

Up to 40°C

Electrical Shock Protection
Class

Class III (per IEC standard)

Massager Type

Roller

Auto Shut-Off

**After 10 minutes of continuous
operation**

Charging Time

Up to 4 hours

Battery Life

**Up to 75 minutes (depending
on mode and load)**

USB Type-C Cable Length

1 m

Device Dimensions

130 × 16.8 × 17.3 cm

Weight

1.2 kg

PACKAGE CONTENTS:


Massager - 1 pc.

USB Type-C cable - 1 pc.

Instruction manual - 1 pc.

DESCRIPTION OF DESIGN AND CONTROLS

EN

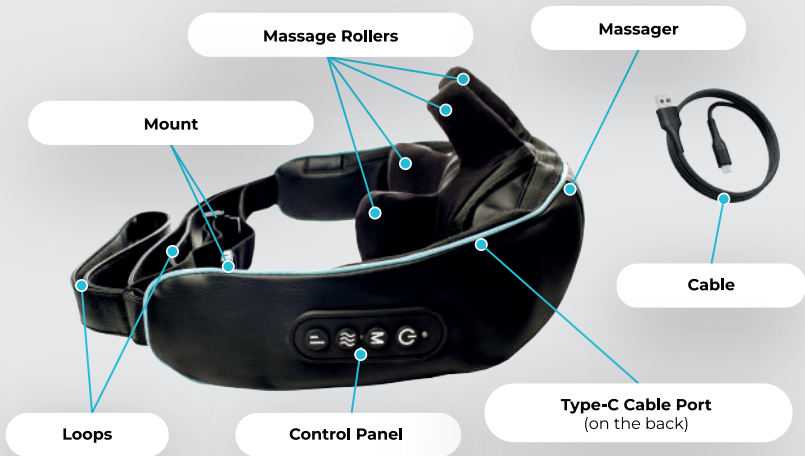DESCRIPTION OF CONTROL BUTTONS

01 Operation and Charging Indicator

Indicates the massager's operation status and charging state

02 Power On/Off Button

Press and hold for 2 seconds to turn the massager on or off

03 Massage Mode Button

Switches between continuous and intermittent modes

04 Heating Button

Turns the heating function on or off

05 Intensity Button

Switches between intensity levels

DESCRIPTION OF CONTROL BUTTONS

DURING CHARGING:

Flashing red

Battery is charging

Solid green for 10 seconds, then turns off

Battery is fully charged

Default Settings Upon Power-On

Message Mode: Continuous

Heating: On

Intensity: Low

DURING MESSAGE:

Solid red

Device is operating, heating function is active

Solid green

Device is operating, heating function is off

Flashing red

Low battery, device will shut off soon

USAGE

- ◆ Press and hold the **«Power On/Off»** button for 2 seconds to start the massager

If needed, press:

- ◆ **«Massage Mode»** — to change the massage type
- ◆ **«Heating»** — to turn the heating function on or off
- ◆ **«Intensity»** — to select the desired intensity level

- ◆ Session duration is — **10 minutes.**
After the session, the device can be manually turned off or will automatically shut off.

IMPORTANT

If any components are missing or the device is damaged, contact the seller or manufacturer for instructions regarding repair or replacement.

SAFETY PRECAUTIONS

EN

GENERAL SAFETY REQUIREMENTS:

- ◆ If you have any medical conditions, consult a doctor before using the device.
- ◆ Always follow the manufacturer's recommendations and safety guidelines to avoid injury or harm during use.
- ◆ Carefully read the user manual before operating the massager.
- ◆ Do not massage the front of the neck (trachea and carotid artery areas).

RECOMMENDATIONS FOR PROPER MASSAGE:

- ◆ Warm up your muscles with light stretching or a warm shower before the massage
- ◆ Massage only the sides and back of the neck, avoiding direct pressure on the spine
- ◆ Use the massager with moderate pressure, especially on sensitive areas
- ◆ Massage different areas of the neck for even muscle relaxation
- ◆ Monitor your sensations: Stop the session immediately if you experience pain, dizziness, or numbness in the limbs
- ◆ Limit each session to 10 minutes, with a maximum of 20 minutes per day
- ◆ Allow your body a few minutes to rest after the massage

INDICATIONS AND CONTRAINDICATIONS

INDICATIONS:

- ◆ Tension, pain, and fatigue in the neck area
- ◆ Recovery after prolonged sitting or static work

CONTRAINDICATIONS:

- ⊗ Arterial hypertension (stages II-III)
- ⊗ Cardiovascular diseases
- ⊗ Blood circulation disorders
- ⊗ Oncological conditions
- ⊗ Pregnancy without a doctor's consultation
- ⊗ Presence of medical implants in the neck area
- ⊗ Skin conditions: acne, eczema, or psoriasis in the treatment area
- ⊗ State of alcohol or drug intoxication
- ⊗ Elevated body temperature or acute inflammatory processes

SAFETY PRECAUTIONS

SPECIAL SAFETY PRECAUTIONS:

- ◆ Use only original accessories and components
- ◆ Do not use the massager on open wounds, tumors, or fresh scars (less than 6 months old)
- ◆ Avoid use by children without a doctor's recommendation
- ◆ Do not use the massager near water, in humid environments, or close to heat sources
- ◆ Store the massager in a dry, safe place, out of reach of children
- ◆ Do not use the massager when extremely fatigued or drowsy

WARNINGS

CAUTION! IF YOU EXPERIENCE PAIN, DIZZINESS, OR ANY DISCOMFORT DURING USE OF THE MASSAGER, STOP THE SESSION IMMEDIATELY AND CONSULT A DOCTOR.

CAUTION! AVOID USING THE MASSAGER ON SENSITIVE AREAS OF THE FACE, AROUND THE EYES, MOUTH, AND NOSE.

POSSIBLE ISSUES AND SOLUTIONS

THE PRODUCT DOES NOT TURN ON

Possible Actions

Check if the device is charged

Press and hold the "Power On" button for 2 seconds to turn on the massager

THE DEVICE DOES NOT CHARGE, AND THE INDICATOR FLASHES RED

Possible Actions

Check if there is power in the outlet

If the USB Type-C cable is damaged, replace it with a functioning one

THE PRODUCT SHUTS OFF DURING USE

Possible Actions

The product automatically shuts off after 10 minutes

Overheating protection has been triggered. Wait 30 minutes, then try turning the device on again

The battery is low, indicated by a flashing red light. Charge the battery

If the issue persists, contact the seller. Do not use the massager until it has been inspected by professionals

POSSIBLE ISSUES AND SOLUTIONS

EN

A SLIGHT NOISE OR CREAKING IS HEARD DURING OPERATION

Possible Actions

The massager's mechanism naturally produces a slight noise during operation. Noise may also occur due to friction between the mechanism and the fabric. This is a normal part of the device's operation and does not require any action

If your issue is not listed above, please contact us at info@zabobon.ua, including photos or videos documenting the problem. We will contact you to help resolve the issue

Read the instruction manual

Recycle damaged or defective products in accordance with local and applicable regulations

Approved in accordance with relevant directives

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